



Clinical Microbiology

Evaluation of a selective chromogenic medium for detecting vancomycin-resistant enterococci



Renata Oliveira Soares*, Adriana Medianeira Rossato, Gustavo Enck Sambrano, Neidimar Cezar Corrêa Tolfo, Juliana Caierão, Thiago Galvão da Silva Paim, Pedro Alves d'Azevedo

Federal University of Health Sciences of Porto Alegre, Porto Alegre, RS, Brazil

ARTICLE INFO

Article history:

Received 7 June 2016

Accepted 1 March 2017

Available online 8 June 2017

Associate Editor: Agnes Figueiredo

ABSTRACT

Rapid identification of vancomycin-resistant enterococci (VRE) can assist in choosing the appropriate treatment and preventing VRE spread. The performance of chromID™ VRE agar was evaluated using 184 clinical isolates of *Enterococcus* spp. and reference strains. The test had a sensitivity of 95.52% but a low specificity of 30%.

© 2017 Sociedade Brasileira de Microbiologia. Published by Elsevier Editora Ltda. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Vancomycin-resistant enterococci (VRE) are among the major agents of healthcare-associated infections and are considered a public health problem. Rapid VRE identification can assist in choosing the appropriate treatment and preventing VRE spread.^{1,2} The aim of this study was to evaluate the performance of a selective chromogenic medium for the detection and differentiation of vancomycin-susceptible and -resistant *Enterococcus faecium* and *Enterococcus faecalis*.

Vancomycin-susceptible enterococci (VSE) isolated from cases of infection ($n=50$) and VRE clinical isolates ($n=134$), including those collected from surveillance rectal swab cultures ($n=62$) and from cases of infection ($n=72$), were evaluated (Table 1). The following reference strains were also included: VRE strain (*E. faecium*, $n=1$) and VSE strains (*E. hirae*, $n=1$; *E. gallinarum*, $n=2$; *E. faecium*, $n=1$; and *E. faecalis*, $n=4$).

All isolates were previously identified by phenotypic methods (hydrolysis of esculin in the presence of bile, production of pyrrolidonyl arylamidase, growth in broth containing 6.5%

NaCl, and negative catalase test evidenced by the absence of effervescence).³ Polymerase chain reaction was also used to confirm the presence of the genus *Enterococcus* and distinguish the species according to methods previously described by Ke et al.⁵ and Karyama et al.⁴ All isolates were obtained from the culture collection of the Gram-positive Cocci Laboratory – UFCSPA and stored in skim milk (Difco™) at -20°C . Vancomycin minimum inhibitory concentration was determined by broth microdilution according to CLSI guidelines (2015)⁶ and by Etest® according to the manufacturer's instructions. The chromID™ VRE (bioMérieux, Brazil S/A) assays were performed in two steps. First, the isolates that were previously stored in skim milk were grown in bile-esculin agar to confirm the presence of enterococci and check culture purity. Second, pure samples were grown in trypticase soy agar for 24 h, followed by single-colony growth in chromID™ VRE agar at 37°C . After 24 h, plates with any growth were considered VRE positive. A negative result was defined as a 48-h incubation period

* Corresponding author.

E-mail: reosoares@gmail.com (R.O. Soares).


<http://dx.doi.org/10.1016/j.bjm.2017.03.005>

1517-8382/© 2017 Sociedade Brasileira de Microbiologia. Published by Elsevier Editora Ltda. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Table 1 – Clinical isolates used to evaluated the performance of the chromID™ VRE.

Species	Susceptibility to vancomycin ^a		Origin	
	S	R	Infection	Surveillance culture
E. faecalis	43	47	47	43
E. faecium	7	87	75	19
Total	50	134	122	62

^a S, susceptible; R, resistant.

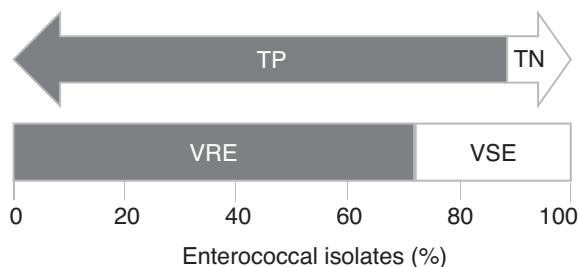

TP = Test positive, TN = Test negative

Fig. 1 – Schematic representation of the results obtained using chromID™ VRE.

without any bacterial growth. According to the manufacturer's instructions, chromID™ VRE agar allows the identification of species based on the detection of enzyme activity. Therefore, *E. faecium* was stained purple and *E. faecalis* was stained blue-green.

Of the *E. faecalis* VRE tested ($n=47$), 41 were stained blue-green, 5 were stained gray, and 1 isolate did not grow. All *E. faecium* VRE tested ($n=87$) were stained purple. Among VSE isolates ($n=50$), 15 did not grow (10 *E. faecalis* and 5 *E. faecium*) and 35 (33 *E. faecalis* and 2 *E. faecium*) showed some growth at the edges of the plates in the corresponding color of each species, which may suggest false-positive results. All VSE-reference strains tested ($n=8$) did not show any visible growth in the chromogenic medium. Fig. 1 shows a schematic representation of the results obtained with chromID™ VRE agar.

The chromID™ VRE agar had a sensitivity of 87.23% and 100% for detecting *E. faecalis* VRE and *E. faecium* VRE, respectively, and a combined sensitivity and specificity of 95.52% and 30.00%, respectively, for detecting VRE. No difference was observed in the specificity and sensitivity at 24 and 48 h. The positive predictive value (corresponding to the percentage of VRE that tested positive in chromID™ VRE) was 78.53% (95% confidence intervals; C.I. = 72.22–84.83%) and the negative predictive value was 71.43% (95% C.I. = 52.10–90.75%) (Table 2). Regarding sensitivity, similar results have been obtained in previous studies evaluating the performance of chromID™ VRE. ^{1,7–14} In relation to the specificity, previous studies have obtained values higher than 95%, ^{1,7,8,11} in contrast to the low specificity observed in this study.

Colonies with non-discriminatory staining (gray, dark, or colorless) have been reported in some studies. ^{10,15} We also observed VRE colonies with a grayish shade, which may lead to false-negative results.

Among all VRE isolates, 100% of the *E. faecium* grew within 24 h. The same was observed in previous studies. ^{11,14} According to Grabsch et al., ⁸ 24-h identification of VRE allows earlier confirmation of colonization by these strains, facilitates infection control, and helps to avoid the spread of microorganisms. However, one *E. faecalis* VRE isolate did not grow even after 48 h, which may suggest that some strains can exhibit a different behavior and/or require more time to grow in the medium.

Delmas et al. ¹ compared growth before and after an enrichment step in bile-esculin agar supplemented with vancomycin in order to select only VRE strains. The enrichment step improved the performance of chromID™ VRE at 24 h of incubation. Other studies have shown that strains incubated overnight in an enrichment broth containing vancomycin as a first step followed by the use of chromID™ VRE resulted in improved specificity or sensitivity. ^{1,10,11,14,16} However, this method is only useful for fecal specimens due to the large number of different microorganisms that can be present in these samples.

Most studies evaluating chromID™ VRE performance have used only fecal specimens (stool samples and rectal swabs) or only resistant strains. In our study, we evaluated well-characterized vancomycin-resistant and -susceptible isolates in order to observe the possible occurrence of false-positive results, because incorrectly prescribed antibiotics have a negative clinical impact. Despite the data presented here, a possible limitation of this study is that fecal samples were not included,

Table 2 – Evaluation of chromID™ VRE in detecting true positive vancomycin-resistant enterococci.

ChromID™ VRE agar	Gold standard (VRE) ^a	Gold standard (VSE) ^a	Total
Test positive	128	35	163
Test negative	6	15	21
Total	134	50	184
Sensitivity		95.52% (95% C.I. = 92.02–99.02%)	
Specificity		30.00% (95% C.I. = 17.29–42.70%)	
Positive predictive value		78.53% (95% C.I. = 72.22–84.83%)	
Negative predictive value		71.43% (95% C.I. = 52.10–90.75%)	
Accuracy		77.72% (95% C.I. = 71.70–83.73%)	

^a Vancomycin resistance was determined using both Etest® and CLSI guidelines; C.I., confidence intervals.

because the density of microorganisms in a clinical specimen may affect the correct diagnosis.

In conclusion, the chromID™ VRE agar is a rapid and useful tool for the screening and identification of VRE, with a good sensitivity of about 96.00%. However, because specificity (30.00%) was limited by false positive VRE mainly, we recommend further VRE identification by conventional tests to avoid misinterpretation.

Conflicts of interest

The authors declare no conflicts of interest.

Acknowledgments

This study was supported in part by the Coordenação de Aperfeiçoamento de Pessoal de Nível Superior and Universidade Federal de Ciências da Saúde de Porto Alegre. We also thank the bioMérieux, Brazil S/A for the supply of the chromID™ VRE agar plates.

REFERENCES

1. Delmas J, Robin F, Schweitzer C, Lesens O, Bonnet R. Evaluation of a new chromogenic medium, chromID VRE, for detection of vancomycin-resistant enterococci in stool samples and rectal swabs. *J Clin Microbiol.* 2007;45:2731–2733.
2. Orsi GB, Ciorba V. Vancomycin resistant enterococci healthcare associated infections. *Ann Ig.* 2013;25:485–492.
3. Teixeira LM, Carvalho MG, Shewmaker PL, et al. Enterococcus. In: Versalovic J, Carroll KC, Funke G, et al., eds. *Manual of Clinical Microbiology*. Washington, DC: ASM Press; 2011:350–364.
4. Karyama R, Mitsuhashi R, Chow JW, Clewell DB, Kumon H. Simple and reliable multiplex PCR assay for surveillance isolates of vancomycin-resistant enterococci. *J Clin Microbiol.* 2000;38:3092–3095. <http://www.ncbi.nlm.nih.gov/pmc/articles/PMC87194/>.
5. Ke D, Picard FJ, Martineau F, et al. Development of a PCR assay for rapid detection of enterococci. *J Clin Microbiol.* 1999;37:3497–3503. <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC85677/>.
6. CLSI – Clinical and Laboratory Standards Institute. Performance Standards for Antimicrobial Susceptibility Testing. In: *Twenty-Fourth Informational Supplement (M100-S24)*, vol. 34. 2015:124–139.
7. Cuzon G, Naas T, Fortineau N, Nordmann P. Novel chromogenic medium for detection of vancomycin-resistant *Enterococcus faecium* and *Enterococcus faecalis*. *J Clin Microbiol.* 2008;46:2442–2444.
8. Grabsch EA, Ghaly-Derias S, Gao W, Howden BP. Comparative study of selective chromogenic (chromID VRE) and bile esculin agars for isolation and identification of vanB-containing vancomycin-resistant enterococci from feces and rectal swabs. *J Clin Microbiol.* 2008;46:4034–4036.
9. Asir K, Wilkinson K, Perry JD, Reed RH, Gould FK. Evaluation of chromogenic media for the isolation of vancomycin-resistant enterococci from stool samples. *Lett Appl Microbiol.* 2009;48:230–233.
10. Peltroche-Llacahuanga H, Top J, Weber-Heynemann J, Lütticken R, Haase G. Comparison of two chromogenic media for selective isolation of vancomycin-resistant enterococci from stool specimens. *J Clin Microbiol.* 2009;47:4113–4116.
11. Lee SY, Park KG, Lee GD, Park JJ, Park YJ. Comparison of Seeplex VRE detection kit with ChromID VRE agar for detection of vancomycin-resistant enterococci in rectal swab specimens. *Ann Clin Lab Sci.* 2010;40:163–166.
12. Klare I, Fleige C, Geringer U, Witte W, Werner G. Performance of three chromogenic VRE screening agars, two Etest® vancomycin protocols, and different microdilution methods in detecting vanB genotype *Enterococcus faecium* with varying vancomycin MICs. *Diagn Microbiol Infect Dis.* 2012;74:171–176.
13. Gouliouris T, Blane B, Brodrick HJ, et al. Comparison of two chromogenic media for the detection of vancomycin-resistant enterococcal carriage by nursing home residents. *Diagn Microbiol Infect Dis.* 2016;85:409–412.
14. Seo JY, Kim PW, Lee JH, et al. Evaluation of PCR-based screening for vancomycin-resistant enterococci compared with a chromogenic agar-based culture method. *J Med Microbiol.* 2011;60(Pt 7):945–949.
15. Ledebuer NA, Tibbets RJ, Dunne WM. A new chromogenic agar medium, chromID VRE, to screen for vancomycin-resistant *Enterococcus faecium* and *Enterococcus faecalis*. *Diagn Microbiol Infect Dis.* 2007;59:477–479.
16. Suwantarar N, Roberts A, Prestridge J, et al. Comparison of five chromogenic media for recovery of vancomycin-resistant enterococci from fecal samples. *J Clin Microbiol.* 2014;52:4039–4042.